

TONTARRA



**Repair and
Exchange Service
available**

FLEXIBLE ENDOSCOPES

Flexible Uretero-Renoscope | Flexible Cystoscope


Our flexible fiberscopes for Urology provide an excellent image quality paired with an ergonomic handle design which results in less fatigue during the intervention.

Furthermore the fiberscopes have a slightly offset deflection tip, thus when using the endoscopes with a laser fiber it is made sure that the laser fiber does not point directly onto the endoscopes' tube.

Such design in combination with a robust mechanical construction and resistant material, ensures a greater durability and extended life-time of the scope.

TONTARRA
made by
EMOS TECHNOLOGY GmbH

TONTARRA


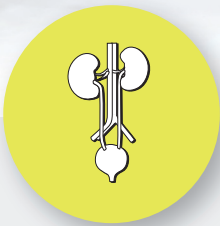


Automatic Valve (Duckbill)
for Uretero-Renoscope;
1 box = 10 pcs.
Art.No. 499-196-00

FLEXIBLE URETERO-RENSCOPE

Insertion Tube	Ø 2.5 mm / Ø 7.5 Fr.
Working channel	Ø 1.2 mm
Working length	680 mm
Field of view	88°
Angulation up/down	275° / 275°
Art.No.	EMTO/FURS-2812R

TONTARRA

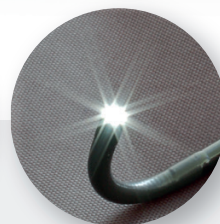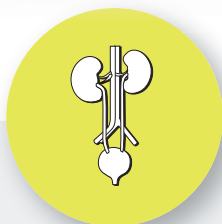

FLEXIBLE CYSTOSCOPE

Insertion Tube	Ø 5.5 mm
----------------	----------

Working channel	Ø 2.3 mm
-----------------	----------

Working length	380 mm
----------------	--------

Field of view	100°
---------------	------

Angulation up/down	210° / 120°
--------------------	-------------

Art.No.	EMTO/FCS-5523R
---------	----------------

- ▶ Available with or without suction valve
(Please specify with order)

TONTARRA
made by
EMOS TECHNOLOGY GmbH



BRILLIANT IMAGE QUALITY

High-quality components „made in Germany“
High depth of field and homogeneous illumination



RESISTANT SHAFT AND GUIDING PART

Extended lifetime due to innovative and robust material and mechanical construction method



SLIGHTLY OFFSET DEFLECTION TIP

Helps avoiding damages of the sheath when using a laser fiber.



FOCUSABLE

Optimal image setting through the ergonomic focusing ring



UNIVERSAL LIGHT GUIDE / CABLE CONNECTION

Compatible with adapters from different manufacturers



STERILIZABLE

with EtO-Gas and Steris®

MADE IN GERMANY



TONTARRA Medizintechnik GmbH

Daimlerstr. 15 · 78573 Wurmlingen/Germany

Tel: +49(0) 74 61/9 65 76 -0

Fax: +49(0) 74 61/9 65 76 -26

info@tontarra.de

www.tontarra.de

All rights are reserved. Technical alterations are possible. This folder may be used for no other purpose than offering, buying and selling of our products. No part may be copied or reproduced in any form. In the case of misuse, we retain the rights to recall our catalogs and pricelists and to take legal actions.

999-405-01E · Rev. 02/24.10.2022